

Titan



Auto



Manual



My Favorite



Position
Adjustment




Dynamic
Squeeze



Setting



Current Massage Program

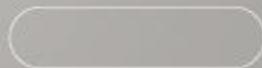



Therapy



DEMO



THAI STRETCH



BODY
CIRCULATION



RELAX



MORNING
CARE



MIDDAY CARE




NIGHT CARE



DEEP TISSUE



Current Massage Program

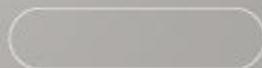



Target



EXTENSION



WOMEN'S
HEALTH



SHIATSU



MEN'S HEALTH



OFFICE RELIEF



SPORTS
RECOVERY



FOOT CARE



REJUVENATE



Current Massage Program





Focus




NECK
MASSAGE



NECK AND
SHOULDER
MASSAGE



BACK
MASSAGE

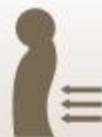

LUMBAR
MASSAGE



WAIST HIP
MASSAGE



LEGS
MASSAGE



UPPER BODY



LOWER BODY



Current Massage Program






Body Detection


Massage chair will automatically detect the body and adjust the seat position. Please keep your head resting on the headrest and sit in the deepest position.

The massage system is going to detect the body.....

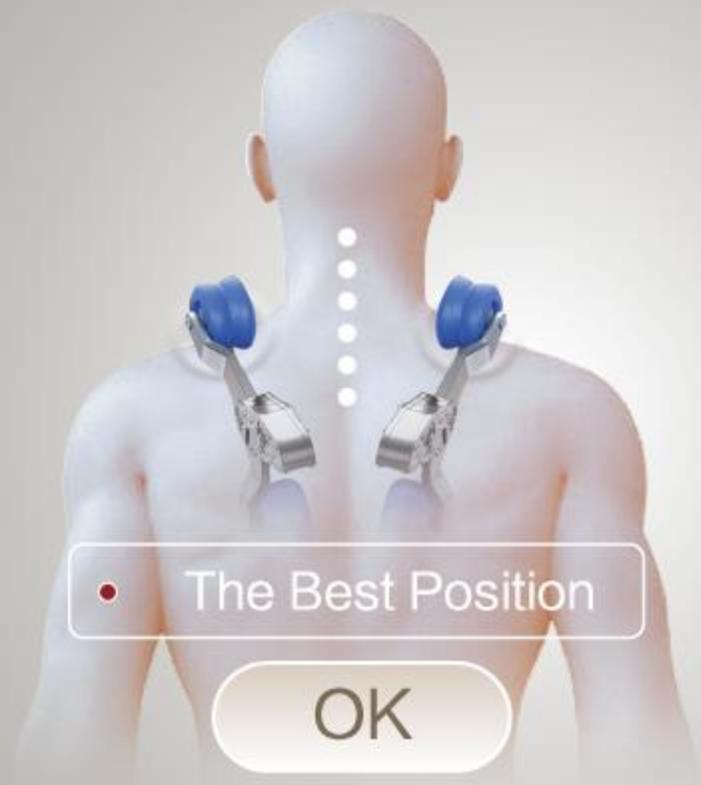

- The Best Position

OK

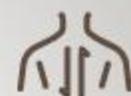 **Shoulder Adjust**

Please press up and down to
adjust the roller position...



Auto

Manual

Dynamic Squeeze

Position Adjustment

★ My Favorite




Core Strength +

Air Intensity +

3D

Air Intensity -

Core Strength -

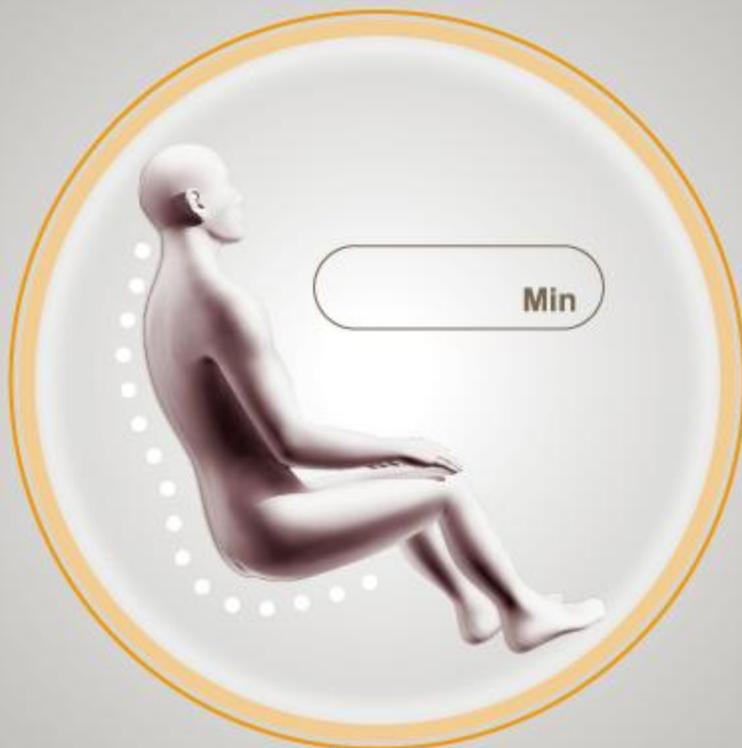

Massage Speed +

Massage Time +

Massage Time -

Massage Speed -



Roller Position



Heating



Current Massage Program





Manual



Auto

Manual


Dynamic Squeeze

Position Adjustment

★ My Favorite



Roller Position Adjustment

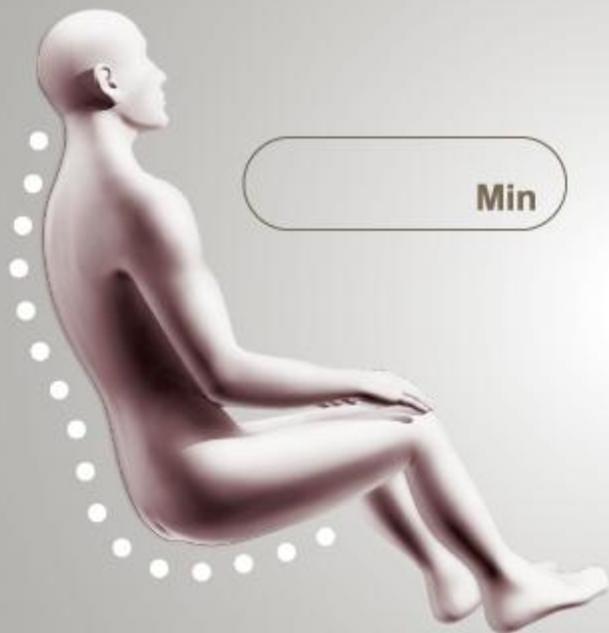

Roller Width

Area Selection



Kneading Kneading & Tapping Knock Tapping Shiatsu Combination




Roller Position



Heating



Current Massage Program

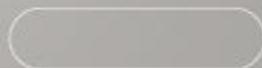

Dynamic
Squeeze



Auto

Manual

Dynamic
Squeeze

Position
Adjustment

★ My Favorite

1. 2. 3. 4. 5.

— Air Intensity

Min

Air Intensity +


Full Body Shoulder Arm Leg & Feet



Current Massage Program



★ My Favorite

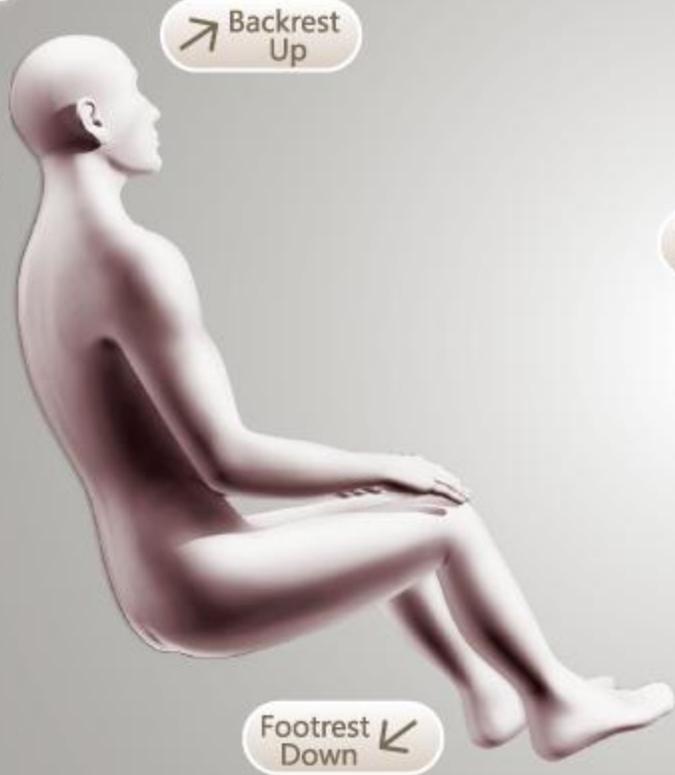

↗ Backrest Up

↙ Backrest Down

↖ Foot Shrink

↘ Foot Extend

↙ Footrest Down

↗ Footrest Up

Zero Gravity



Watching TV Reading



⚠ When adjusting position, children sitting in front or behind the massage chair, please note do not allow pets to approach.

Current Massage Program





Current Massage Program





Function Setting

Language

Others



Leg Kneading




Foot Roller



Heating



Current Massage Program

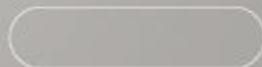



Function Setting

Language

Others


English

한국어

简体中文



Current Massage Program

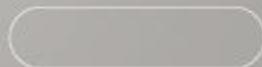



Function Setting


Language

Others

Voice Commands



Dormancy Time

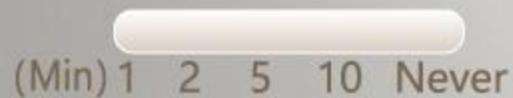

Voice Control



Backlight



Key Tone



Bluetooth

Sleep Mode

Current Massage Program



Voice control function

- Go to " SETTING" and tap on the "VOICE CONTROL" icon button to activate, then say " Hi,Titan/ Hey,Titan ", you will get a reply "What can I do for you?"
- The commands are shown in table below.

	Voice Commands		Voice Commands		Voice Commands		Voice Commands
1	HI,TITAN / HEY,TITAN	10	NIGHT CARE MODE	19	REJUVENATE MODE	28	INCREASE STRENGTH
2	PAUSE MASSAGE	11	DEEP TISSUE MODE	20	NECK MASSAGE MODE	29	DECREASE STRENGTH
3	STOP MASSAGE	12	EXTENSION MODE	21	NECK AND SHOULDER MASSAGE MODE	30	INCREASE AIR
4	DEMO MODE	13	WOMEN'S HEALTH MODE	22	BACK MASSAGE MODE	31	DECREASE AIR
5	THAI STRETCH MASSAGE MODE	14	SHIATSU MODE	23	LUMBAR MASSAGE MODE		
6	BODY CIRCULATION MODE	15	MEN'S HEALTH MODE	24	WAIST & HIP MASSAGE MODE		
7	RELAX MASSAGE MODE	16	OFFICE RELIEF MODE	25	LEGS MASSAGE MODE		
8	MORNING CARE MODE	17	SPORTS RECOVERY MODE	26	UPPER BODY MASSAGE MODE		
9	MIDDAY CARE MODE	18	FOOT CARE MODE	27	LOWER BODY MASSAGE MODE		

- The massage chair will execute the commands within 10 seconds.
- If the commands were not said out clearly, the massage chair will not execute the commands.
- If there is no new command after 15 seconds, the voice control system will response "Thank you for using"
- To resume voice control, say " Hi,Titan/Hey,Titan " again.





The massage chair is being reset, please wait...